



### **Amico Group of Companies**

Amico Corporation, founded in 1974, designs, manufactures and markets a broad range of products for virtually every department in the health care facility. Our mission is to provide state-of-the-art products that facilitate the caregiver's role in providing the best outcomes for patients.

These products are manufactured in eight facilities in the U.S. and Canada, which occupy over **700,000** square feet.












Amico Group of Companies - One Point of Contact. One Total Solution.

### **GoLift** Features

- World's lightest patient lift
- Emergency stop
- Emergency lowering
- Battery status indicators
- LED light indicators
- Wired magnetic hand control
- Infection control features
- Manual emergency lowering
- Quick connect carry bar
- Power Traverse option

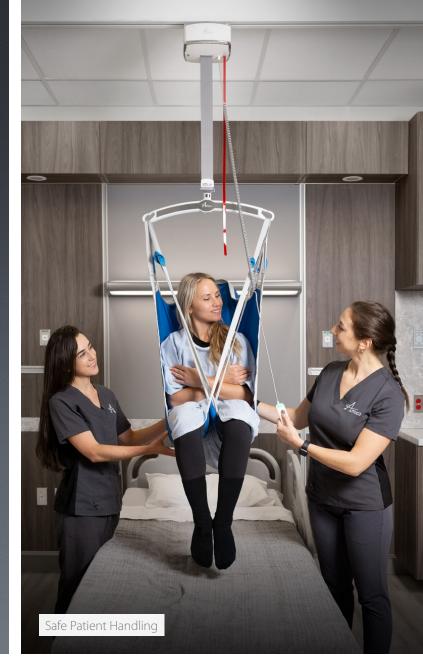










### GoLift Patient Lift Pendant (PLP)



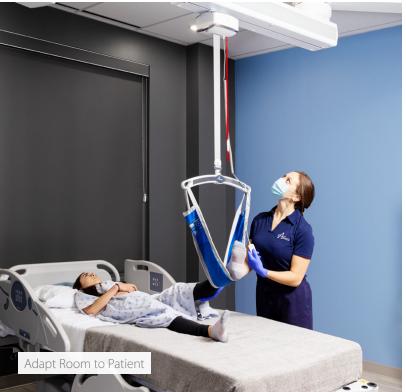






# GoLift Patient Lift Pendant (PLP)

in the OR

. . . .

PLP Outside the Surgical Field

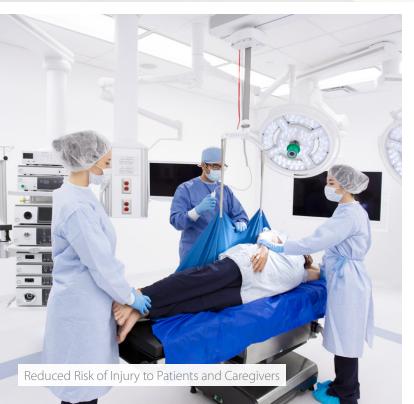



:

i bea

- Base




### GoLift Headwall Integration







### **GoLift** Service & Maintenance

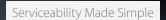





Amico

## Retrofit with GoLift

- Quick connection design
- Adapter trolleys work with all major lift manufacturer tracks
- Easy installation and maintenance
- Compatible with end stop and continuous charging systems
- Instant electrical connection









### **GoLift** 1000

Introducing the GoLift 1000 - a powerful product catered to Bariatric patients. Two compact GoLift 700 systems designed for an increased lifting capacity of 1000 lbs (454 kg).

- Four Point Carry Bar
- Emergency stop
- Emergency lowering
- Battery status indicators
- LED light indicators
- Wired magnetic hand control
- Infection control features
- Manual emergency lowering
- Quick connect carry bar







Four Point Carry Bar



# **GoLift**1000

The GoLift1000 can be converted into two GoLift700<u>s!</u>

GoLift1000 Trolley




### **GoLift** Goes Portable!

The Amico GoLift **Portable450** system is designed for clinicians to transfer a patient lift from room to room. It is the smallest and lightest portable patient lift in the world today.

- Weighs 6 lbs, lifts 450 lbs!
- Wired magnetic hand control
- Controls on lift motor
- Emergency stop and lowering
- Manual lowering
- Reacher Bar
- Long lasting Ni-Mh battery



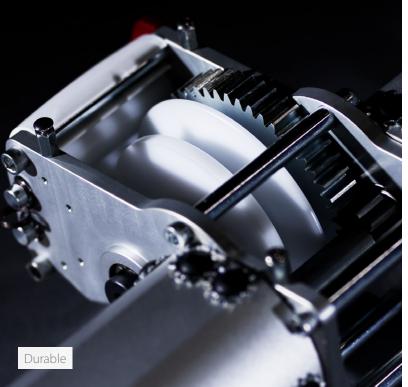







# **GoLift** Portable **450**








## **GoLift** Slings

We offer a wide selection of reusable and disposable slings. Refer to our Sling Brochure for more details.


#### Guidelines for Choosing your Sling Size

				S		М		L		XL
<b>Go</b> Basic		With Head Support	Height: Shoulder: Hip:	29" (73.6 cm) 20" (50.8 cm) 16" (40.6 cm)	Height: Shoulder: Hip:	30" (76.2 cm) 24" (60.9 cm) 20" (50.8 cm)	Height: Shoulder: Hip:	31" (78.7 cm) 28" (71.1 cm) 24" (60.9 cm)	Height: Shoulder: Hip:	32" (81.2 cm) 32" (81.2 cm) 28" (71.1 cm)
		No Head Support	Height: Shoulder: Hip:	21" (53.3 cm) 20" (50.8 cm) 16" (40.6 cm)	Height: Shoulder: Hip:	22" (55.8 cm) 24" (60.9 cm) 20" (50.8 cm)	Height: Shoulder: Hip:	23" (58.4 cm) 28" (71.1 cm) 24" (60.9 cm)	Height: Shoulder: Hip:	24" (31 cm) 32" (81.2 cm) 28" (71.1 cm)
Golong			42" x 78" (106.6 cm x 198.1 cm)			56" x 78" (142.2 cm x 198.1 cm)				
GHygiene			Chest:	32" (81.2 cm)	Chest:	36" (91.4 cm)	Chest:	40" (101.6 cm)	Chest:	44" (111.7 cm)
<b>GOCOMFORT</b>			Height: Shoulder: Hip:	23" (58.4 cm) 24" (60.9 cm) 20" (50.8 cm)	Height: Shoulder: Hip:	24" (60.9 cm) 26" (66 cm) 22" (55.8 cm)	Height: Shoulder: Hip:	25" (63.5 cm) 28" (71.1 cm) 24" (60.9 cm)	Height: Shoulder: Hip:	26" (66 cm) 30" (76.2 cm) 26" (66 cm)
GOActive	Classic General Constants Classic General Constants		Chest:	40" (101.6 cm)	Chest:	44" (111.7 cm)	Chest:	48" (121.9 cm)	Chest:	52" (132 cm)
<b>Go</b> Reposition						20" x 75" (51 c	m x 190.5 cr	n)		
GoLimb						15" x 7" (38.1	cm x 17.8 cm	1)		
GoTurn	X					Width: 22'	" (55.8 cm)			
<b>GO</b> Amputee	-		Height: Hip:	28" (71.1 cm) 15" (38.1 cm)	Height: Hip:	33" (83.8 cm) 18" (45.7 cm)	Height: Hip:	38" (96.5 cm) 22" (55.9 cm)	Height: Hip:	43" (109.2 cm) 26" (66 cm)
Gowalk			Hip:	15" (38.1 cm)	Hip:	18" (45.7 cm)	Hip:	22" (55.9 cm)	Hip:	26" (66 cm)

### Anatomy of the GoLift

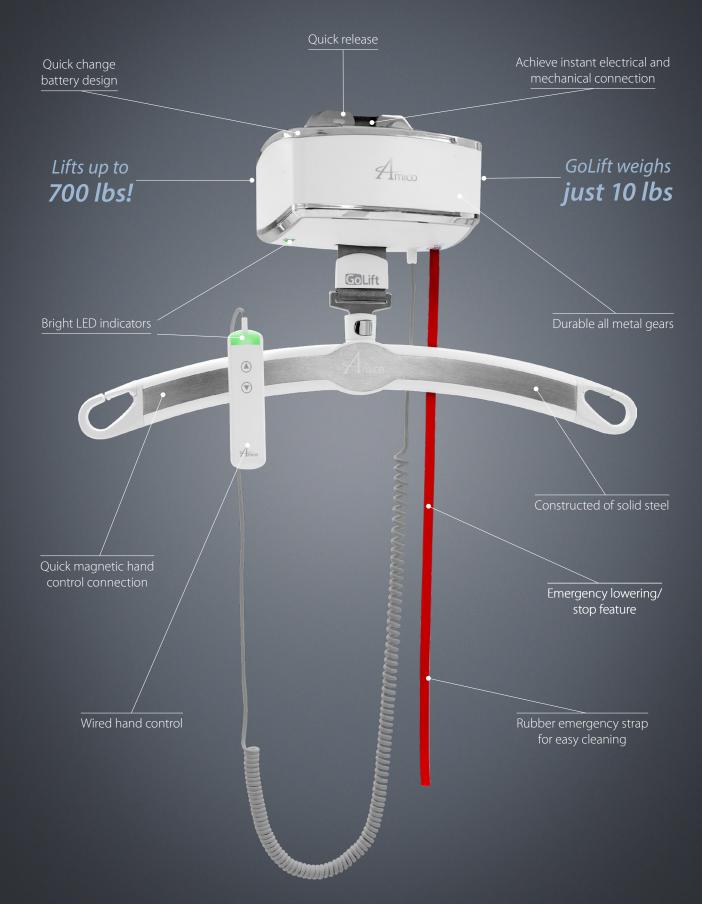

### Anatomy of the GoLift Portable450

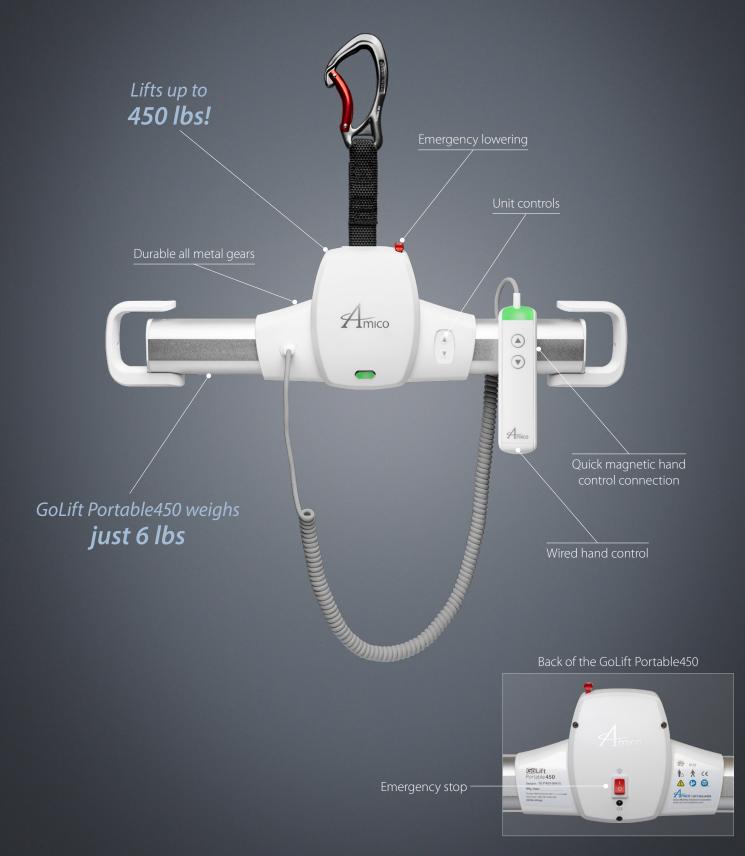

### GoTrack Systems


### Configurations

- XY Gantry (H-Track)
- Straight Track
- Patient Lift Pendant
- Turntable (Manual or Power)
- Free-Standing Track
- Wall-Post Track
- Power Traverse Gantry


### Features

Straight Track

- Light weight and high lifting
- Quick and easy
- Hidden end stops for safety
- Continuous charge or end stop charge

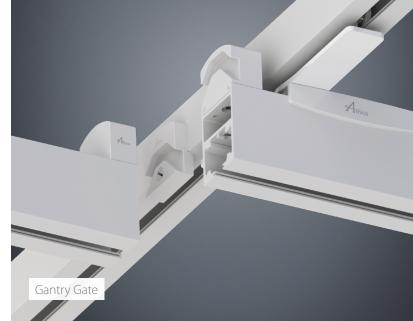



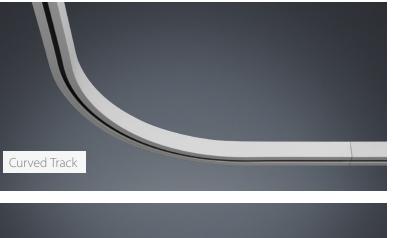



### **GoLift** Accessories



One Hook Carry Bar



Two Hook Carry Bar



Four Point Carry Bar



Fixed One Hook Carry Bai



Fixed Two Hook Carry Bar



Power Traverse



Curtain Jumping Trolley



GoLift400 Trolley



GoLift700 Trolley





### **GoLift** Specifications



#### **Maximum Lifting Speed**

No load:2 inches/secondGoLift400:1 inch/secondGoLift700:1 inch/secondGoLift1000:1 inch/second

#### **Unit Weight**

 GoLift400:
 10 lbs (4.5 kg)

 GoLift700:
 10 lbs (4.5 kg)

 GoLift1000:
 20 lbs (9 kg) (Tandem Lifts)

#### Safe Working Load (SWL)

### GoLift400:400 lbs (181 kg)GoLift700:700 lbs (318 kg)GoLift1000:1000 lbs (454 kg)

#### Dimensions

GoLift400 / GoLift700:

٠	Length:	7.75"	(196.8 mm)	
	VA /: -L+L-	7701	(1000)	

• vviatri.	/./ 3 (190.0 11111)
Lloight	4.25''(1070 mm)

- Height: 4.25" (107.9 mm)
  Strap Length: 84" (2133.6 mm)
- Strap Length. 04 (2155.01111)

#### GoLift1000:

- Length: 15.5" (393.7 mm)
- Width: 7.75" (196.8 mm)
- Height: 4.25" (107.9 mm)
- Strap Length: 84" (2133.6 mm)

#### Number of Lifts per Charge (Duty: 10/90)

25% of Strap at Midrange<sup>1</sup>:

- 375 with 185 lbs (84 kg)
- 170 with 400 lbs (181 kg)
- 102 with 700 lbs (318 kg)

• 140 with 1000 lbs (454 kg) Charging Time: 2-4 hours *Calculated using 5Ah Packs* 

#### Lift Case

Flame Retardant ABS

#### **Batteries**

High Capacity, Nickel Metal Hydride (NiMH)Standard:2x 14.4V (3.3Ah)Optional:2x 14.4V (5Ah)

#### **Battery Charger:**

 Input:
 100-240V AC

 Output:
 36V DC, 1.0A, 40W

#### Hand Control

Protection Class: IPX4

**Maximum Sound Level** 

#### **Expected** Service Life

10 years or 22,500 cycles

#### Safety

Emergency Stop Emergency Lowering Device Upper Limit Detection Lower Limit Detection Slack Tape Sensor Free Fall Brake (over speed governor) Low Battery and Dead Battery Alarms Soft Start and Stop Overload Protection Emergency Manual Lowering

#### Approvals

Certified to:

Tested to:

CAN/CSA-C22.2 No 60601-1:08, UL 60601-1:08 ISO 10535-06, CE







#### **Maximum Lifting Speed**

No load: 150 lbs/68 kg: 450 lbs/204 kg:

#### 1.7 inches/second 1.3 inches/second 1.0 inch/second

#### Unit Weight

6 lbs (2.72 kg)

#### Safe Working Load (SWL)

450 lbs (204 kg)

#### Dimensions

- Length: 19.4" (492.7 mm)
- Height: 5.2" (132 mm)
- Depth of Lift: 3.7" (94 mm)
- Strap Length: 74" (1879.6 mm)

#### Number of Lifts per Charge (Duty: 10/90)

- 25% of Strap at Midrange:
- 78 with 150 lbs (68 kg)
- 46 with 450 lbs (204 kg)

Charging Time: 2-4 Hours

#### Lift Case

Flame Retardant ABS IP23 Lift - protected from fingers and water spray less 60 degrees from vertical

#### Batteries

High Capacity, Nickel Metal Hydride (NiMH) 1x 24V (1.6Ah)

#### **Battery Charger**

 Input:
 100-240V AC

 Output:
 36V DC, 1.0A, 40W

#### Hand Control

IPX4 Hand Control - protects from splashing water, no matter the direction

#### Expected Service Life

22,500 cycles or 10 years. It is recommended that the lift be replaced after this time. The expected life of slings, batteries fuses, straps and cord is dependent upon proper maintenance.

#### Safety

Emergency Stop Emergency Lowering Device: Independent Circuit Manual Lowering Upper Limit Detection Lower Limit Detection Slack Tape Sensor Free Fall Brake (Over Speed Governor) Low Battery and Dead Battery Alarms Soft Start and Stop Overload Protection

#### Approvals

CAN/CSA-C22.2 No 60601-1:08 CAN/CSA-Z10535.1:15





For free survey / quotes and brochures Please call us: 843-547-2257



InMyHomeAccess@gmail.com www.InMyHomeAccess.com